

ELECTROENCEPHALOGRAPHY AND REACTION TIME-Recording

Visual Reaction Times

Casandra Ortega

Jasmin Villagomez

Purpose

Our bodies respond to various sensory inputs. This procedure was used to examine the time it took for my lab partner and I to respond to a visual reaction. While we know that there is a reaction time, we wanted to test it ourselves and see how long it took for her and I to react.

Procedures

1. To get things started:

Before you turn anything on, be sure the IWX/214 unit is plugged in, and that the IWX/214 unit is connected to the laptop by USB cable.

Be sure that the EM-100Event Marker is fully connected to the Channel 3 socket in front of the IWX/214 unit (Fig. 5-4)

Once everything is connected, FIRST turn on the laptop and allow it to fully boot up before you turn on the IWX/214 unit. Once the Iworx unit is on, the red indicator light on the Iworx unit should light up and you may hear the USB chime from the laptop if the laptop does not default to mute (many are set to default to mute)

2. Open the Labscribe 3 Program by clicking on the LabScribe3iconon the desktop. As soon as the program opens, you should see a window pop-up that says "Hardware found IWX 214:2008-1-24," click "OK." If you see a window pop-up with procedures for installing new hardware, you may have turned on the IWX/214 unit before the laptop fully booted up. Simply turn off the IWX/214 unit, close the Labscribe 3 Program and start over again by first turning the IWX/214 back on and then reopening the Labscribe 3 Program as described above

3. In the second from the top row (the row that says "File Edit View Tools Settings AdvancedExternal DevicesHelp"), click on the "Settings" tab. A little more than halfway down the drop-down window should be a tab called "Human Nerve." Click on that tab and that should lead you to the second tab called "Auditory-VisualReflexes," click on that tab. Close the pdf file that opens automatically, you don't need it.

4. Pair up with a lab partner and arrange yourselves according to this recording setup: The student subject should sit in a chair facing the laptop computer with their hand in position so that they can press the “Enter” key as quickly as possible. The lab partner holding the EM-100 Event Marker Should stand out of sight of the subject. They need to be able to quietly press and release the button of the Event Marker once the test begins

5. Type the subject’s name and “Visual” in the Mark box that is to the right of the Mark button just above the data recording. Click the red “Record” button then click the Mark button, this will put a vertical line in your recording and the words in the Mark box at the bottom of the vertical line. Leave the cursor mouse over the Mark box (not Mark button)

6. Each time the lab partner quietly clicks the Event Marker button, the green line coming in from the right side of the computer screen will jump up then back down. As soon as the subject sees the green line jump up, they need to press the “Enter” key as quickly as possible. When this happens a small “Enter Mark Text ” window will pop up, ignore it. However, the subject will have to click on the “Enter” key again to clear it. Do not worry, the Mark line will still be in the right place. The subject will have to hit the “Enter” key twice for every trial: once to leave the Mark line and a second time to clear the “Enter Mark Text” pop up window.

7. Repeat this for ten trials but the lab partner should be sure to click the Event Marker button at irregular intervals (not less than 5 seconds apart, but not more than 10 seconds apart). After ten trials click the “Stop” button (it’s the same button that turns from the red “Rec” button to the black “Stop” button once it is recording).


8. Multiple lab partners can use the same Iworx file. Just click the red “Rec” button again and let a good 20 seconds of flat line go by before you click the black “Stop” button. Then repeat steps 5-7 above with the new subject’s name. If a third student will be using the same file, just click the red “Rec” button again and let a good 20 seconds of flat line go by before you click the black “Stop” button. Then repeat steps 5-7 with the third student’s name.

9. Once all lab partners have made their Visual Cues recording, go back to the first of ten trials for each student. Move the red cursor lines (there are two) by left clicking on the red cursor line and while holding down the left touchpad button, sliding your finger on the mouse pad to move one red cursor line to the left base of the green rectangle, then releasing the left touchpad button. Move The other red cursor line to the black Mark line.

10. Look at the top right of the screen where it says “ $T_2 - T_1$ ” = _____ msec. Record that number and repeat this step for all ten trials and for each of the students using your laptop. Calculate the average for the ten trials. Report These numbers to the lab instructor who will compile the class numbers.

Results

Jas vs. Cass

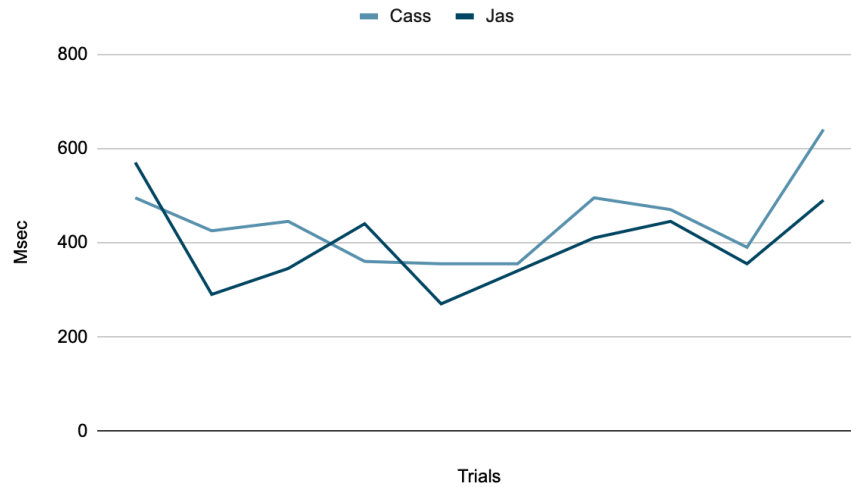

Cass	495	425	445	360	355	355	495	470	390	640
Jas	570	290	345	440	270	340	410	445	355	490

Discussion

While ours are not exactly the same, the results were still very similar. The average for Cass was 407.5 msec and the average for Jasmin was 395.5 msec. It was very close in range showing that our results are normal. One second is not a long time, it is clear that our reaction time is very quick according to our results.

Conclusion

In conclusion, we examined how fast our response was according to our visual input to the green box shown on the screen. Our average was very similar showing that our results are normal.